

SECRO-FILL

Body Filler


2×50ml pre-filled syringe

Made in Germany





Product Type	Biphasic Cross-Linked Hyaluronic Acid Filler
Volume	2 × 50ml (Total 100ml, 2 pre-filled syringes)
HA Concentration	26mg/ml
Core Application	Buttock Augmentation & Contouring: Enhances buttock volume for a fuller, lifted shape; refines contour and improves symmetry, creating a harmonious body curve.
Suggested Injection Depth	Subcutaneous to deep subcutaneous tissue

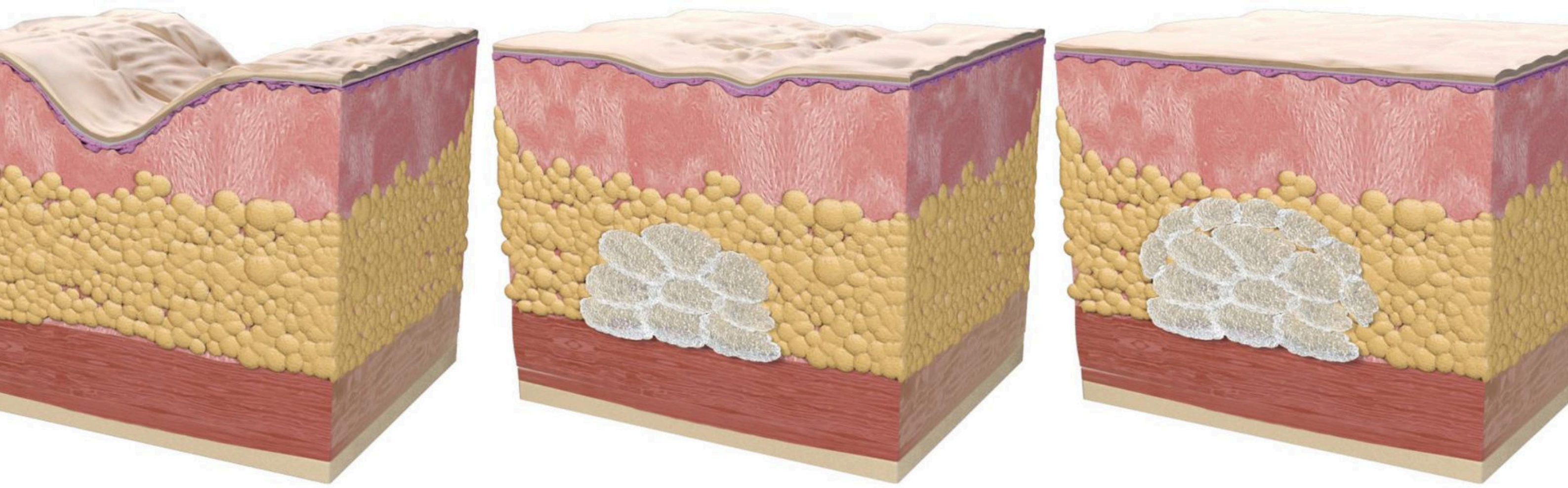

Filling skin depressions
BUTT ENLARGEMENT



What are the advantages of SECRO-FILL BODY fillers?

- **2-syringe pre-filled packaging:**
Ready-to-use, no complicated transfer or dosage measurement steps
- Boosts clinical efficiency and reduces contamination risk during preparation
- Flexible dosage adjustment for personalized buttock enhancement needs
- **26mg/ml** high concentration for strong & stable support; long-lasting effect up to 12-18 months; fully biocompatible with minimal downtime



BIPHASIC GEL STRUCTURE ADVANTAGE:

Combines a firm cross-linked HA particle core for durable volume retention with a smooth non-cross-linked HA matrix for natural tissue integration, balancing long-term support and soft, natural post-treatment texture.

FEEDBACKS:



CONTRAINDICATIONS & PRECAUTIONS

Contraindications:

- Patients with hypersensitivity to hyaluronic acid or any ingredients in the product
- Patients with active infections, inflammation, or skin lesions at the injection site
- Pregnant or breastfeeding women
- Patients with bleeding disorders or taking anticoagulant medications

Precautions:

- For professional use only by certified aesthetic physicians.
- Do not inject intravascularly to avoid vascular occlusion.
- Post-treatment care: Avoid strenuous exercise, saunas, or facial/body massages for 24 hours; apply cold compresses to reduce swelling if needed.



SECROMED

📍 HAMBURG, GERMANY

☎ +491 516 708 3711

✉ Info@secromed.com

🌐 www.secromed.com

📱 @secromed.com

