



SECROMED

HYAC-LIFT (16% CHAC)

(1 × 2ml pre-filled syringe)





80 mg high molecular CHAC HA + 80 mg low molecular CHAC HA



SECROMED
PROFESSIONAL USE

HYAC-LIFT (16% CHAC)

Pre-filled syringe. Net content 2 mL
Microbiologically tested
Hyaluronic acid 160 mg/mL





CHAC technology is able to modify and utilize the natural properties of hyaluronic acid, making it the best carrier for delivering essential nutrients to the skin. The mixture under goes pressure and shear deformation simultaneously, which creates a molecular complex.

The additives are mostly small organic molecules that stabilize the hyaluronic acid, thereby increasing the bioavailability of the additives and providing long-lasting effects.

HYAC Injectable Boosters are made with the highest quality hyaluronic acid, amino acids and vitamins, etc. These ingredients are carefully selected to treat various skin problems and effectively fight aging and make the skin younger.

They are 100% pure active complexes without compromising safety.

What is CHAC Technology?

HYAC
PROFESSIONAL USE

What is HYAC-LIFT (16% CHAC)?

HYAC-LIFT (16% CHAC) is a very effective anti-aging product.

Its formula contains different molecular weights of hyaluronic acid obtained through a bio-fermentation process:

80 mg high molecular weight CHAC HA + 80 mg low molecular weight CHAC HA.

No chemical modification is required. Hyaluronic acid acts as a humectant, helping to promote hydration in all skin types, and is especially beneficial for dry and dehydrated skin.

It helps rejuvenate the outer surface layer of the skin by increasing the production of elastin, fibronectin, aminoglucan and collagen, making the skin look and feel softer, smoother and more hydrated. By increasing elastin and fibronectin, HA functions in the extracellular matrix to create physiological conditions for the proliferation, migration and organization of the dermal cell components.

HYAC
PROFESSIONAL USE

Why are there different molecular weights of hyaluronic acid in the same solution?

A 2016 study answered this question:


The permeability of hyaluronic acid of different molecular weights in human skin was probed by Raman spectroscopy.

The results were very clear:

Material with molecular weights up to 300 kDa has good penetration,

materials between 50-300 kDa are ideal for anti-aging, while materials with lower molecular weights (80-10 kDa) can penetrate deeper and are ideal for eliminating cellulite and retaining moisture.

On the other hand, substances with molecular weights above 1000 kDa stay on the surface of the epidermis, while the body's natural 2000 kDa hyaluronic acid is completely unable to penetrate the epidermis.

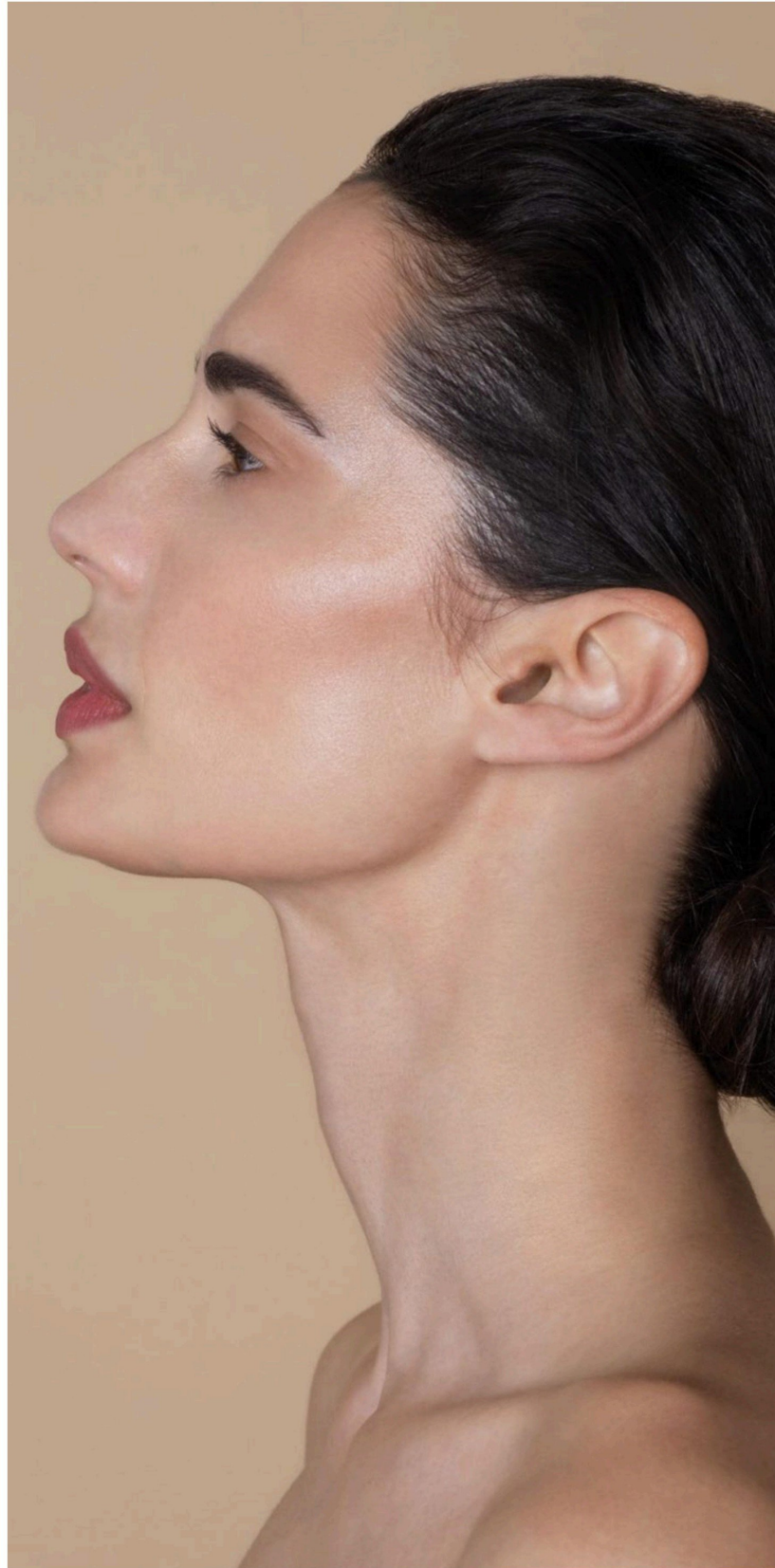

Main ingredients of HYAC-LIFT (16% CHAC)

80 mg high molecular CHAC HA

+

80 mg low molecular CHAC HA



Sodium
chloride

CHAC
Sodium
hyaluronate

Sodium
phosphate

Disodium
phosphate

HYAC
PROFESSIONAL USE

The Effects of HYAC-LIFT (16% CHAC):

Moisturizing: Deeply hydrates, nourishes the skin, regulates the water-oil balance of the skin, and improves skin problems caused by dehydration.

Lifting and firming: Quickly locks deep moisture in the skin and promotes the production of subcutaneous collagen and elastin. Lifts loose skin to make the facial contour appear fuller, smoother and more beautiful.

Restore elasticity: Promote the regeneration of subcutaneous collagen and elastin, restore skin elasticity, and restore a youthful look.

Anti-aging and anti-wrinkle: Improve cell vitality, smooth fine lines and wrinkles. Regulate the cell microenvironment and improve skin quality.



What are the advantages of HYAC-LIFT (16% CHAC)?

- The product has excellent fluidity and ductility, and the skin and the product can blend naturally, which makes the skin very natural and not stiff after injection.
- The injection course is convenient and fast, and you can immediately resume normal work and life after surgery.
- Double skin rejuvenation effect: Help the skin restore its youthful state of tenderness, smoothness and firmness, moisturize and lock in moisture, tighten the skin, and remove wrinkles.
- High concentration: The effect exceeds that of ordinary hyaluronic acid, the effect is long-lasting and the process is safe.



HYAC
PROFESSIONAL USE

Five-point confirmation:

- Protruding cheekbone: at least 2 cm from the outer corner of the eye.
- Bottom of the nose: draw a line connecting the nostril and the tragus, draw a vertical line from the pupil, and position the injection point at the intersection of the two lines.
- Tragus: 1 cm in front of the bottom of the tragus.
- Chin: Draw a vertical line in the center of the chin, draw a vertical line from the top of the vertical line to the lower third, and move 1.5 cm from the intersection point to the oral commissure.
- Mandibular angle: 1 cm upward from the mandibular angle.

HYAC
PROFESSIONAL USE





HYAC
PROFESSIONAL USE

Applicable parts:

Face, neck, shoulders, arms, hands, knees

Applicable people:

This product is suitable for people of all skin types, especially for those over 25 years old with fine lines and slight sagging.

Number of treatments:

2 sessions every 1 month

Precautions:

- Before treatment, please keep your skin clean, and any cosmetics and skin care products need to be clean.
- It is strongly recommended that patients do not take anti-inflammatory drugs (such as aspirin, etc.) at least 3-4 days before treatment, as it may aggravate bleeding and swelling at the injection site.



SECROMED
Made in Germany

www.secromed.com